

इंटरनेट

मानक

Disclosure to Promote the Right To Information

Whereas the Parliament of India has set out to provide a practical regime of right to information for citizens to secure access to information under the control of public authorities, in order to promote transparency and accountability in the working of every public authority, and whereas the attached publication of the Bureau of Indian Standards is of particular interest to the public, particularly disadvantaged communities and those engaged in the pursuit of education and knowledge, the attached public safety standard is made available to promote the timely dissemination of this information in an accurate manner to the public.

“जानने का अधिकार, जीने का अधिकार”

Mazdoor Kisan Shakti Sangathan

“The Right to Information, The Right to Live”

“पुराने को छोड़ नये के तरफ”

Jawaharlal Nehru

“Step Out From the Old to the New”

IS 3025-48 (1994): Methods of sampling and test (physical and chemical) for water and wastewater, Part 48: Mercury [CHD 32: Environmental Protection and Waste Management]



“ज्ञान से एक नये भारत का निर्माण”

Satyanarayan Gangaram Pitroda

“Invent a New India Using Knowledge”



“ज्ञान एक ऐसा खजाना है जो कभी चुराया नहीं जा सकता है”

Bhartrhari—Nitiśatakam

“Knowledge is such a treasure which cannot be stolen”

BLANK PAGE



IS 3025 (Part 48) : 1994

(Reaffirmed 1999)

REAFFIRMED

२००३

भारतीय मानक

जल और अपशिष्ट जल के नमूने लेने (भौतिक एवं रसायन)
तथा परीक्षण की पद्धति

भाग 48 पारा

(पहला पुनरीक्षण)

Indian Standard

METHODS OF SAMPLING AND TEST
(PHYSICAL AND CHEMICAL) FOR WATER
AND WASTEWATER

PART 48 MERCURY

(*First Revision*)

Second Reprint JANUARY 2002

UDC 628.1.032 : 628.3 : 543.3 [546.49]

© BIS 1994

BUREAU OF INDIAN STANDARDS
MANAK BHAVAN, 9 BAHADUR SHAH ZAFAR MARG
NEW DELHI 110002

March 1994

Price Group 3

AMENDMENT NO. 1 OCTOBER 2000
TO
IS 3025 (PART 48) : 1994 METHODS OF SAMPLING
AND TEST (PHYSICAL AND CHEMICAL) FOR WATER
AND WASTEWATER
PART 48 MERCURY

(First Revision)

(Second Cover page, Foreword, line 12) — Substitute the following for the existing matter:

'This standard supersedes 13 of IS 2488 (Part 2) : 1968 Methods of sampling and test for industrial effluents Part 2.'

(Page 3, clause 6.4.7, line 10) — Substitute 'ammonium hydroxide' for 'minimum hydroxide'.

(Page 3, clause 6.4.7, line 23) — Insert the words 'completely. Add 25 ml of chloroform, shake well' after the words 'dithizone precipitates'.

(Page 4, clause 6.5.1, line 19) — Insert the words 'with 10 ml chloroform solution' after the word 'extraction'.

(Page 4, clause 6.5.1, line 36) — Insert the word 'chloroform layer' after the words 'mercury dithizone'.

(CHD 12)

FOREWORD

This Indian Standard (Part 48) (First Revision) was adopted by the Bureau of Indian Standards, after the draft finalized by the Environmental Protection Sectional Committee had been approved by the Chemical Division Council.

Organic and inorganic mercury salts are very toxic and their presence in the environment, especially in water, should be monitored. Therefore, it is desirable to ascertain the limit of mercury present in water and wastewater. In the preparation of this standard, considerable assistance has been derived from Standard Methods for Examination of Water and Wastewater, 16th Edition, 1985.

The technical committee responsible for formulation of IS 2488 and IS 3025 : 1964 decided to revise all the parts covered under them and publish as separate parts.

This standard supersedes 7 of IS 2488 (Part 5) : 1976 'Methods of sampling and test for industrial effluents, Part 5.

The composition of the technical committee responsible for the formulation of this Indian Standard is given in Annex A.

In reporting the result of a test or analysis made in accordance with this standard, if the final value, observed or calculated, is to be rounded off, it shall be done in accordance with IS 2 : 1960 'Rules for rounding off numerical values (revised)'.

Indian Standard

METHODS OF SAMPLING AND TEST (PHYSICAL AND CHEMICAL) FOR WATER AND WASTEWATER

PART 48 MERCURY

(First Revision)

1 SCOPE

This standard prescribes two methods for the determination of mercury in water and wastewater.

- a) Cold vapour atomic absorption spectrophotometric method, and
- b) Calorimetric dithizone method.

2 REFERENCES

The following Indian Standards are necessary adjuncts to this standard:

<i>IS No.</i>	<i>Title</i>
7022 (Part 1) : 1973	Glossary of terms relating to water, sewage and industrial effluents, Part 1
7022 (Part 2) : 1979	Glossary of terms relating to water, sewage and industrial effluents, Part 2

3 TERMINOLOGY

For the purpose of this standard, definitions given in IS 7022 (Part 1) : 1973 and IS 7022 (Part 2) : 1979 shall apply.

4 SAMPLE PRESERVATION

The sampling bottles should be cleaned thoroughly with dilute nitric acid (6 N) prior to final rinsing with water. The water samples should be collected and stored preferably in polypropylene or chemically resistant glass containers. For preservation, the sample should be acidified with concentrated nitric acid (2 ml of AR grade nitric acid in one litre of the sample just to bring down the pH below 2). For dissolved mercury, filter the sample in the field and acidify the filtrate with nitric acid to a pH of 2 or lower.

NOTE — Avoid excess nitric acid. Add 5 ml of 0.1 N iodine solution to avoid losses of volatile organo mercury components during handling and digestion of samples.

5 COLD VAPOUR ATOMIC ABSORPTION SPECTROMETRY

5.1 Principle

The flameless atomic absorption procedure is a physical method based on the absorption of radiation at 253.7 nm by mercury vapour. The mercury is reduced to the elemental state and aerated from solution in a closed system. The mercury vapour passes through a cell positioned in the light path of mercury hollow cathode lamp of an atomic absorption spectrophotometer. Absorbance (peak height) is measured as a function of mercury concentration and recorded.

The cold vapour atomic absorption spectrophotometric method is suitable for all types of samples such as natural waters, potable waters, domestic and industrial wastewaters. Lower detection limit of 0.2 µg/l may be achieved.

5.2 Apparatus

5.2.1 Atomic Absorption Spectrometer (AAS) and Associated Equipment

Instrument settings recommended by the manufacturer shall be followed. Instruments designed specifically for the measurement of mercury using the cold vapour technique may be substituted for the AAS.

5.2.2 Mercury Vapour Generation Assembly

Consists of an absorption cell, peristaltic pump, flow meter, aeration tubing and a drying tube containing magnesium perchlorate.

5.2.3 Mercury Hollow Cathode Lamp

5.2.4 Recorder / Printer / Display Meter

Any multi-range variable recorder that is compatible with the UV detection system is suitable.

5.3 Reagents

5.3.1 Sulphuric Acid — Concentrated.

5.3.2 Nitric Acid — Concentrated.

5.3.3 Stannous Chloride Solution

Dissolve 25 g of stannous chloride (SnCl_2) in water containing 50 ml of concentrated hydrochloric acid and dilute to 250 ml. If a suspension forms, stir reagent continuously during use.

5.3.4 Sodium Chloride — Hydroxylamine sulphate solution.

Dissolve 12 g of sodium chloride and 12 g of hydroxylaminesulphate ($(\text{NH}_2\text{OH})_2\text{H}_2\text{SO}_4$) in distilled water and dilute to 100 ml.

5.3.5 Potassium Permanganate Solution

Dissolve 5 g of potassium permanganate in distilled water and dilute to 100 ml.

5.3.6 Potassium Persulphate Solution

Dissolve 5 g of potassium persulphate in distilled water and dilute to 100 ml.

5.3.7 Stock Mercury Solution

Dissolve 1.354 g of mercuric chloride in about 700 ml of distilled water. Add 10 ml of concentrated nitric acid and make up to 1 000 ml (1 ml = 1 mg Hg).

5.3.8 Standard Mercury Solution

Prepare a series of standard mercury solutions containing 0 to 5 $\mu\text{g/l}$ by appropriate dilution of stock mercury solution (5.3.7) with water

containing 10 ml of concentrated nitric acid per litre. Prepare standards daily.

NOTE — Use mercury free distilled water for the preparation of reagents and standards.

5.4 Procedure

5.4.1 Instrument Operation

Follow the procedure of the manufacturer's operating manual. Connect the mercury vapour generating assembly as shown in Fig. 1.

5.4.2 Standardization

Transfer 100 ml of each of the 1.0, 2.0 and 5.0 $\mu\text{g/l}$ standard mercury solution and a blank of 100 ml water to 300 ml BOD bottles. Add 5 ml of concentrated sulphuric acid and 2.5 ml of concentrated nitric acid to each bottle. Add 15 ml of potassium permanganate solution to each bottle and let stand for at least 15 minutes. Add 8 ml of potassium persulphate ($\text{K}_2\text{S}_2\text{O}_8$) solution to each bottle and heat for 2 hours in a water bath at 95°C. Cool and add 6 ml of sodium chloride-hydroxylamine sulphate solution to reduce the excess permanganate. After decolourization add 5 ml of stannous chloride solution and attach the bottle immediately to the aeration apparatus forming a closed system. As mercury is volatilised and carried into the absorption cell, absorbance will increase to a maximum within a few seconds. As soon as recorder returns approximately to the base line, remove stopper holding the aeration frit from the reaction bottle and replace with a bottle containing distilled water. Flush the system for a few seconds and run the next standard in

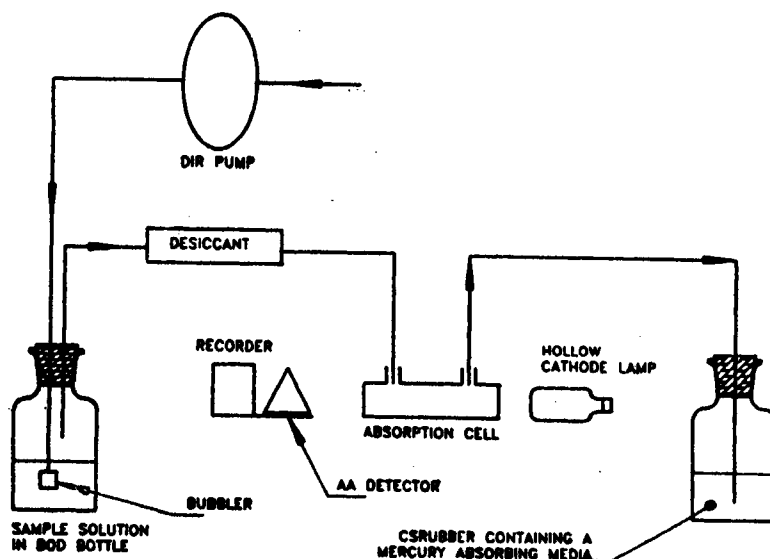


FIG. 1. SCHEMATIC ARRANGEMENT OF EQUIPMENT FOR MEASUREMENT OF MERCURY BY COLD VAPOUR ATOMIC ABSORPTION TECHNIQUE

the same manner. Construct a standard calibration curve by plotting absorbance (peak height) versus mercury concentration in μg .

5.4.3 Analysis of Sample

Transfer 100 ml of sample or portion diluted to 100 ml containing not more than $5.0 \mu\text{g/l}$ mercury to a 300 ml BOD bottle. Treat as in 5.4.2 samples containing high (sea waters, brines and effluents high in chloride) chlorides require as much as an additional 25 ml potassium permanganate solution. During this step, chloride is converted into free chlorine which is absorbed at 253.7 nm. Hence remove all free chlorine before the mercury is reduced and swept into the cell by using an excess (25 ml) of hydroxylamine sulphate solution.

5.5 Calculation

Determine peak height of sample from recorder chart and read mercury value from standard curve.

6 COLORIMETRIC DITHIZONE METHOD

6.1 Principle

Mercury ions react with dithizone solution in chloroform to form an orange colour which is measured at 490 nm using a spectrophotometer. The colorimetric dithizone method is suitable for estimating higher levels of mercury in potable waters (more than $2 \mu\text{g/l}$) and has a minimum detection limit of $2 \mu\text{g/l}$.

6.2 Interference

Copper, gold, palladium, divalent platinum and silver react with dithizone in acid solution. Copper in the dithizone extract remains in the organic phase while the mercury dissolves in the aqueous phase. The other contaminants usually are not present. The mercury dithizonate being photosensitive should be measured quickly.

6.3 Apparatus

6.3.1 Spectrophotometer

For use at 492 nm, providing a light path of 1 cm or longer.

6.3.2 Separating Funnels

250 and 1 000 ml with PTFE stopcocks.

6.3.3 Glassware

Clean all glassware with potassium dichromate-sulphuric acid cleaning solution.

6.4 Reagents

6.4.1 Mercury Free Distilled Water

Use redistilled or deionised distilled water for preparing all reagents and dilutions.

6.4.2 Stock Mercury Solution

Dissolve 135.4 mg of mercuric chloride (HgCl_2) in about 700 ml of distilled water, add 1.5 ml concentrated nitric acid and make up to 1 000 ml with distilled water ($1.0 \text{ ml} = 100 \mu\text{g Hg}$).

6.4.3 Standard Mercury Solution

Dilute 10.0 ml of stock solution to 1 000 ml with distilled water. Prepare freshly before use ($1.00 \text{ ml} = 1.00 \mu\text{g Hg}$).

6.4.4 Potassium Permanganate Solution

Dissolve 5 g of potassium permanganate in 100 ml of distilled water.

6.4.5 Potassium Persulphate Solution

Dissolve 5 g of potassium persulphate ($\text{K}_2\text{S}_2\text{O}_8$) in 100 ml distilled water.

6.4.6 Hydroxylamine Hydrochloride Solution

Dissolve 50 g of hydroxylamine hydrochloride ($\text{NH}_2\text{OH} \cdot \text{HCl}$) in 100 ml water.

6.4.7 Dithizone Solution

Dithizone is often contaminated with the oxidation product diphenyl thiocarbodizone or with metals. To purify dithizone dissolve 100 mg of dithizone in 50 ml of chloroform in a beaker and filter under slight vacuum using Whatman No. 42 or equivalent filter paper. Wash beaker and filter paper with 5 ml portions of chloroform. Transfer the filtrate to 500 ml separatory funnel. Add 100 ml of 1 percent minimum hydroxide and shake moderately for 1 minute only. Transfer chloroform layer to 250 ml separatory funnel retaining orange-red aqueous layer in a 500 ml funnel. Repeat extraction process keeping each time chloroform layers separately in 250 ml separatory funnels and using 1 percent a minimum hydroxide. After three extractions, transfer aqueous layer to 500 ml funnel. Discard chloroform layer.

Combine all the extracts in a 500 ml separatory funnel. Add 1 : 1 hydrochloric acid in 2 ml portions, mixing each time till dithizone precipitates and orange-red colour appears with three 25 ml portions chloroform. Dilute the combined extract to 1 000 ml with chloroform; $1 000 \text{ ml} = 100 \mu\text{g}$ dithizone. Dilute 60 ml of this solution with chloroform to 1 000 ml, ($1 \text{ ml} = 6 \mu\text{g}$ dithizone).

6.4.8 Sulphuric Acid — 0.25 N.

Dilute 25 ml of 1 N sulphuric acid to 100 ml with distilled water.

6.4.9 Potassium Bromide Solution

Dissolve 40 g of potassium bromide in 100 ml distilled water.

6.4.10 Chloroform

6.4.11 Phosphate-Carbonate Buffer Solution —

Dissolve 150 g of $\text{Na}_2\text{HPO}_4 \cdot 12\text{H}_2\text{O}$ and 38 g of anhydrous potassium carbonate (K_2CO_3) in 1 litre of distilled water. Extract with 10 ml portions of dithizone until the last portion remains blue. Wash with chloroform to remove excess dithizone.

6.4.12 Sodium Sulphate — Anhydrous.

6.5 Procedure

6.5.1 Preparation of Calibration Curve

Pipette 0 (blank), 2.0, 4.0, 6.0, 8.0 and 10.0 ml of mercury standard solutions into separate beakers. To each beaker, add 500 ml of distilled water or that volume chosen for sample, 1 ml of potassium permanganate solution and 10 ml of concentrated sulphuric acid solution. Stir and boil if necessary, add more potassium permanganate until a pink colour persists. After boiling has ceased continuously add 5 ml of potassium persulphate ($\text{K}_2\text{S}_2\text{O}_8$) solution and let cool for 30 minutes. Add a few drops of hydroxylamine hydrochloride ($\text{NH}_2\text{OH} \cdot \text{HCl}$) solution to discharge the pink colour. Cool and transfer each solution to individual 1 litre separating funnels. Add 25 ml of dithizone solution. Shake the separating funnel vigorously and transfer each organic layer to 250 ml separating funnel. Repeat this extraction at least three times making sure that the colour in

the last dithizone extract is an intense blue as that of the original dithizone solution. Wash accumulated dithizone extracts in the 250 ml separating funnel by shaking with 50 ml of 0.25 N sulphuric acid. Transfer washed dithizonate extract to another 250 ml separating funnel. Add 50 ml of 0.25 N sulphuric acid and 10 ml of potassium bromide solution and shake vigorously to transfer mercury dithizonate from organic layer to aqueous layer. Discard dithizone layer. Wash aqueous layer with a small volume of chloroform and discard the chloroform. Transfer 20 ml of phosphate carbonate buffer solution to each separating funnel and 10 ml standard dithizone solution. Shake vigorously and after separation, transfer the mercury dithizone to beakers. Dry contents with anhydrous sodium sulphate. Transfer mercury dithizonate solution to a cuvette and record absorbance at 492 nm. Plot absorbance versus mercury concentration in μg .

6.5.2 Analysis of Samples

Use 500 ml of sample and prepare an absorbance blank consisting of all reagents. When necessary, filter sample through glass wool into the separating funnel after oxidation step. Complete procedure as described in 6.5.1.

6.6 Calculation

6.6.1 Read mercury from calibration curve.

6.6.2 After obtaining the above value, depending upon the volume chosen for analysis, calculate for 1000 ml and report result in $\mu\text{g/l Hg}$.

ANNEX A

(Foreword)

COMMITTEE COMPOSITION

Environmental Protection Sectional Committee, CHD 012

<i>Chairman</i>	<i>Representing</i>
PROF D. K. BISWAS	Central Pollution Control Board, Delhi
<i>Members</i>	
DR K. R. RANGANATHAN (<i>Alternate</i> to Prof D. K. Biswas)	
ADVISER (PHE)	Ministry of Rural Development
ADDL ADVISER (PHE) (<i>Alternate</i>)	
SHRI S. B. C. AGARWALA	Bharat Heavy Electricals Ltd, Hyderabad
SHRI S. BALAGURUNATHAN (<i>Alternate</i> I)	
SHRI A. K. GUPTA (<i>Alternate</i> II)	
DR A. L. AGGARWAL	National Environmental Engineering Research Institute (CSIR), Nagpur
DR T. CHAKRABARTI (<i>Alternate</i>)	
DR A. ALAM	Indian Council of Agricultural Research, New Delhi
SHRI S. C. AHLUWALIA	National Council for Cement & Building Materials, New Delhi
SHRI A. D. AGNIHOTRI (<i>Alternate</i>)	
SHRI R. K. BANERJEE	Shriram Institute for Industrial Research, Delhi
SHRI P. K. MAIR (<i>Alternate</i>)	
SHRI B. BASU	National Thermal Power Corporation Ltd, New Delhi
DR S. MUKHERJEE (<i>Alternate</i>)	
SHRI V. S. BHATNAGAR	Central Scientific Instruments Organization (CSIR), Chandigarh
DR M. S. N. SRINIVAS (<i>Alternate</i>)	
SHRI S. CHAKRAVORTI	Directorate General Factory Advice Service & Labour Institutes, Bombay
DR M. H. FULEKAR (<i>Alternate</i>)	
SHRI S. DAS	Indian Petrochemicals Corporation Ltd, Vadodara
SHRI M. K. PRABHUDESAI (<i>Alternate</i>)	
DR V. S. GUPTA	National Test House, Calcutta
SHRI D. N. P. SINGH (<i>Alternate</i>)	
DR HARISH CHANDRA	Industrial Toxicology Research Centre (CSIR), Lucknow
SHRI B. K. JAIN	The Fertilizer Association of India, New Delhi
DR (Ms) B. SWAMINATHAN (<i>Alternate</i>)	
SHRI G. K. GURJIA	Thermax Ltd, Pune
DR A. K. WAGLE (<i>Alternate</i>)	
SHRI A. LAHIRI	Hindustan Lever Ltd, Bombay
SHRI B. B. DAVE (<i>Alternate</i>)	
DR W. MADHAVAKRISHNA	Central Leather Research Institute (CSIR), Madras
SHRI S. RAJAMANI (<i>Alternate</i>)	
SHRI S. K. MAIRA	Flakt India Ltd, Calcutta
SHRI A. SAHA (<i>Alternate</i>)	
SHRI R. K. MALHOTRA	Indian Oil Corporation Ltd (R & D Centre), Faridabad
SHRI S. K. JAIN (<i>Alternate</i>)	
SHRI A. N. KALE	Municipal Corporation of Greater Bombay
SHRI V. S. MAHAJAN (<i>Alternate</i>)	
DR P. K. MATHUR	Bhabha Atomic Research Centre (IGCAR), Kalpakkam (TN)
DR P. M. MODAK	Indian Institute of Technology, Bombay
PROF H. VEFRAMANI (<i>Alternate</i>)	
SHRI K. P. NYATI	National Productivity Council, New Delhi
SHRI L. PANEERSELVAM (<i>Alternate</i>)	
PROF B. PADMANABHAMURTHY	Jawaharlal Nehru University, New Delhi
DR T. S. PATEL	National Institute of Occupational Health (ICMR), Ahmadabad
SHRI C. V. RAIYANI (<i>Alternate</i>)	
DR V. V. RAO	Dharmshi Morarji Chemical Co Ltd, Bombay
DR M. ATCHAYYA (<i>Alternate</i>)	

(Continued on page 6)

IS 3025 (Part 48) : 1994

(Continued from page 5)

Members

SHRI P. S. RAMANATHAN
SHRI D. N. V. RAO
SHRI R. J. BUCH (Alternate)
DR S. ROUTH
DR J. C. NIJHAWAN (Alternate)
SHRI P. R. SAMADDAR
SHRI P. K. SEN (Alternate)
SHRI S. C. SHARMA
SHRI R. N. GUPTA (Alternate)
SHRI M. P. SINGH
SHRI N. C. TIWARI (Alternate)
SHRI M. SUBBA RAO
DR T. CHANDINI (Alternate)
SHRI R. M. SUNDARAM
SHRI C. KRISHNA RAO (Alternate)
SHRI SURENDER KUMAR
SHRI R. PARTHASARTHY (Alternate)
SUPERINTENDING ENGINEER
EXECUTIVE ENGINEER (Alternate)
SHRI J. M. TULI
SHRI S. N. CHAKRABARTI (Alternate)
DR R. K. SINGH,
Director (Chem)

Representing

Pesticides Association of India, New Delhi
Tata Chemicals Ltd, Bombay
National Test House, Calcutta
Central Mechanical Engg Research Institute (CSIR), Durgapur
India Meteorological Department, New Delhi
Directorate General of Technical Development, New Delhi
Ministry of Environment & Forests
National Malaria Eradication Programme (DGHS), Delhi
Indian Chemical Manufacturers Association, New Delhi
Panchayat Raj Department, Government of Andhra Pradesh,
Hyderabad
Engineers India Ltd, New Delhi
Director General, BIS (Ex-officio Member)

Member Secretary

SHRI T. RANGASAMY
Joint Director (Chem), BIS

Water Environment Subcommittee, CHD 012 : 01

Convener

DR Y. P. KAKAR

Members

SHRI S. B. C. AGARWALA
SHRI A. K. GUPTA (Alternate)
SHRI A. BASU
SHRI A. K. JINDAL (Alternate)
SHRI M. S. DHINGRA
SHRI V. G. K. NAIR (Alternate)
DR E. K. JAYANARAYANAN
SHRI K. K. MITTU (Alternate)
SHRI S. ISLAM
SHRI F. LAL KANSAL
SHRI S. S. SANGHA (Alternate)
SHRI D. D. KUMTA
DR K. C. PATHAK (Alternate)
PROF K. J. NATH
PROF A. K. ADHYA (Alternate)
DR R. NATH
DR S. RATAN (Alternate)
DR S. R. PANDE
DR M. V. NANOTI (Alternate)
DR P. M. PHIRKE
DR S. R. JOSHI (Alternate)
SHRI S. PRAKASH
SHRI S. S. RAMRAKHYANI (Alternate)
SHRI R. V. RAO
SHRI D. K. KAUSHIK (Alternate)
REPRESENTATIVE
REPRESENTATIVE
DR B. SENGUPTA
DR R. C. TRIVEDI (Alternate)
SUPERINTENDING ENGINEER
EXECUTIVE ENGINEER (Alternate)
SHRI S. R. TAMTA
SHRI K. RAJAGOPALAN (Alternate)
DR P. N. VISWANATHAN

Ministry of Environment & Forests

Bharat Heavy Electricals Ltd, Hyderabad

Thermax Ltd, Pune

Shriram Institute for Industrial Research, Delhi

Mohan Meakin Ltd, Mohan Nagar

Central Pulp and Paper Research Institute, Saharanpur
Punjab Pollution Control Board, Patiala

Tata Chemicals Ltd, Bombay

All India Institute of Hygiene & Public Health, Calcutta

Banaras Hindu University, Varanasi

National Environment Engineering Research Institute (CSIR),
Nagpur

National Environment Engineering Research Institute (CSIR),
Nagpur

Delhi Water Supply & Sewage Disposal Undertaking, New Delhi

Central Water Commission, New Delhi

Ministry of Rural Development
U. P. Jal Nigam, Lucknow
Central Pollution Control Board, Delhi

Panchayat Raj Department, Government of Andhra Pradesh,
Hyderabad

Central Ground Water Board, New Delhi

Industrial Toxicology Research Centre (CSIR), Lucknow

Bureau of Indian Standards

BIS is a statutory institution established under the *Bureau of Indian Standards Act, 1986* to promote harmonious development of the activities of standardization, marking and quality certification of goods and attending to connected matters in the country.

Copyright

BIS has the copyright of all its publications. No part of these publications may be reproduced in any form without the prior permission in writing of BIS. This does not preclude the free use, in the course of implementing the standard, of necessary details, such as symbols and sizes, type or grade designations. Enquiries relating to copyright be addressed to the Director (Publications), BIS.

Review of Indian Standards

Amendments are issued to standards as the need arises on the basis of comments. Standards are also reviewed periodically; a standard along with amendments is reaffirmed when such review indicates that no changes are needed; if the review indicates that changes are needed, it is taken up for revision. Users of Indian Standards should ascertain that they are in possession of the latest amendments or edition by referring to the latest issue of 'BIS Catalogue' and 'Standards: Monthly Additions'.

This Indian Standard has been developed from Doc : No. CHD 012 (0171

Amendments Issued Since Publication

Amend No.	Date of Issue	Text Affected

BUREAU OF INDIAN STANDARDS

Headquarters :

Manak Bhavan, 9 Bahadur Shah Zafar Marg, New Delhi 110 002
Telephones : 323 01 31, 323 33 75, 323 94 02

Telegrams : Manaksanstha
(Common to all offices)

Regional Offices :

	Telephone
Central : Manak Bhavan, 9 Bahadur Shah Zafar Marg NEW DELHI 110 002	{ 323 76 17 323 38 41
Eastern : 1/14 C.I.T. Scheme VII M, V. I. P. Road, Kankurgachi KOLKATA 700 054	{ 337 84 99, 337 85 61 337 86 26, 337 91 20
Northern : SCO 335-336, Sector 34-A, CHANDIGARH 160 022	{ 60 38 43 60 20 25
Southern : C.I.T. Campus, IV Cross Road, CHENNAI 600 113	{ 254 12 16, 254 14 42 254 25 19, 254 13 15
Western : Manakalaya, E9 MIDC, Marol, Andheri (East) MUMBAI 400 093	{ 832 92 95, 832 78 58 832 78 91, 832 78 92
Branches : AHMEDABAD. BANGALORE. BHOPAL. BHUBANESHWAR. COIMBATORE. FARIDABAD. GHAZIABAD. GUWAHATI. HYDERABAD. JAIPUR. KANPUR. LUCKNOW. NAGPUR. NALAGARH. PATNA. PUNE. RAJKOT. THIRUVANANTHAPURAM.	

Reprography Unit, BIS, New Delhi, India